

**ottobock.**

# Agilium Reactive

For unicompartmental osteoarthritis of the knee.  
Less pain. More life.



Quality for life

# Agilium Reactive

## For knee osteoarthritis

Osteoarthritis is the most common joint disorder worldwide. In light of our aging population, the number of those affected is increasing. Women and people who are overweight are diagnosed more frequently than other groups.

As an orthotist, you are familiar with the symptoms experienced by people with knee osteoarthritis: the most basic everyday activities can cause knee pain. On the other hand, these activities can feel great when pain is relieved.



### Did you know?

An osteoarthritis patient takes an average of around 1,200 tablets a year to manage pain. But this can lead to damage to the stomach, bowel and liver. An orthosis from the Agilium line offers a good alternative.<sup>1</sup> It's worthwhile for anyone with knee osteoarthritis to test the effectiveness of the orthoses themselves.

The Agilium Reactive is a knee orthosis that uses the traditional three-point fitting principle for medial compartment knee osteoarthritis. It is optimised for patients, fitting providers and prescribers who are aiming to achieve a comprehensive effect on the knee. Its innovative fastening technology in the upper calf area provides comfort while sitting without preventing stable positioning when standing.

<sup>1</sup> Schmalz et al. 2006; Schmalz et al. 2011; Fantini-Pagani et al. 2014; Stinus et al. 2015 (not yet published)



### For a high level of comfort

The Agilium Reactive is lightweight. This is due to the material used and a slim and flat design. The innovative surface technology of the padding enables comfort when worn yet ensures good protection against slipping.

### Sit more comfortably at the press of a button

Users often have the problem that the orthosis creates uncomfortable pressure when they are sitting down. The calf muscle changes its shape when the user sits down and the upper calf strap suddenly becomes too tight, when it was perfectly comfortable whilst standing. The Agilium Reactive offers a unique solution to this. Thanks to the inbuilt Boa fastening system, the strap can be loosened at the press of a button. If stable support is required again, the user can increase the support with a few turns of the Boa fastening.


### Adapts to the shape of the body

Thanks to the auto-adaptive design of the straps, the Agilium Reactive can be adjusted to fit various leg shapes. This makes it possible to offer the orthosis in only one universal size.



### Sizing chart for 50K324 Agilium Reactive

Order examples: 50K324=L-1  
50K324=R-1

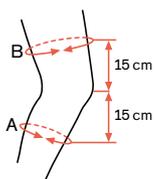

Side	Circumference (cm)	
L/R	A	30–60
	B	39–80



Agilium Freestep 2.0



Agilium Patella

### Agilium line

The orthoses in our Agilium line are designed specifically to target the symptoms of knee osteoarthritis. Each works in a different way to address the various characteristics of knee osteoarthritis. At the same time, we placed special focus on their comfort and suitability for daily use.

